




Peninsula Laboratories International, Inc.

305 Old County Road, San Carlos, CA 94070
Tel: (650) 801-6090
Fax: (650) 595-4071



Cat. No. S-1420 Amylin (human) (for human serum and plasma samples)

SPECIFICATIONS: Extraction-free kit for human serum/plasma
STANDARD: 1 µg H-7905 H-Lys-Cys-Asn-Thr-Ala-Thr-Cys-Ala-Thr-Gln-Arg-Leu-Ala-Asn-Phe-Leu-Val-His-Ser-Ser-Asn-Asn-Phe-Gly-Ala-Ile-Leu-Ser-Ser-Thr-Asn-Val-Gly-Ser-Asn-Thr-Tyr-NH₂, (Disulfide bond)
PROTOCOL: V (Ab1hr.Std2hr.BtON)
RANGE: 0-100 ng/ml
AVERAGE IC₅₀: 1.7 ng/ml
SHELF-LIFE: One year from production date

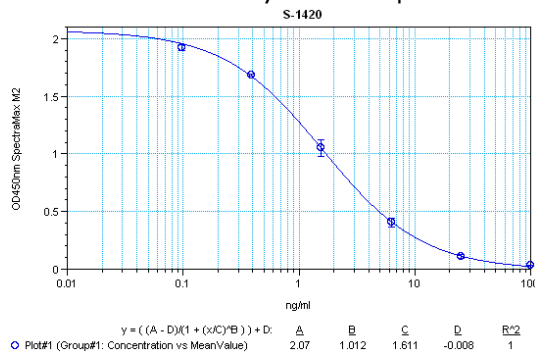

Suggested Preparation of Standards		
	ng/ml	Range: 0.10 to 100 ng/ml
Stock	1000	
S1	100.00	Add 100 µl stock + 900 µl diluent
S2	25.00	Add 200 µl S1 + 600 µl diluent
S3	6.25	Add 200 µl S2 + 600 µl diluent
S4	1.56	Add 200 µl S3 + 600 µl diluent
S5	0.39	Add 200 µl S4 + 600 µl diluent
S6	0.10	Add 200 µl S5 + 600 µl diluent
S0	0.00	500 µl diluent

CROSS REACTION:

PEPTIDE	%
Amylin (free acid) (human)	100
Amylin (human)	100
Amylin (rat)	100
Amylin (feline)	100
α-CGRP (rat)	75
α-CGRP (human)	1
Calcitonin (human)	0
NPY (porcine)	0
Pancreastatin (porcine)	0
VIP (human, bovine, rat, porcine)	0
Insulin (human)	0
Glucagon (1-29) (human, bovine, porcine)	0
Somatostatin-14	0
Pancreatic Polypeptide (human)	0
Proinsulin (human)	0

Cross-Reactivity data is taken from Radioimmunoassay applications.